

Keeler Loupe Frame

INSTRUCTIONS FOR USE

CONTENTS

1. INDICATIONS FOR USE	2
1.1 WARNINGS AND CAUTIONS	2
1.2 CONTRAINDICATION	3
2. CLEANING AND DISINFECTION INSTRUCTIONS	4
3. SETTING UP AND USING THE LOUPE FRAME.....	4
4. TECHNICAL SPECIFICATIONS.....	6
5. SPARE PARTS AND ACCESSORIES	6
6. WARRANTY	6
7. PACKAGING AND DISPOSAL INFORMATION.....	7

	Consult instructions for use		General warning sign
	Date of manufacture		Conformité Européene
	Manufacturer's name and address		Temperature limit
	Country of manufacture		Authorised representative in the European Community
	Waste Electrical and Electronic Equipment (WEEE) recycling		Catalogue number
	This way up		Batch code
	Keep dry		Atmospheric pressure limitation
	Fragile		Humidity limitation
	Do not use if package is damaged		Translation

The Keeler Loupe Frame is designed and built in conformity with Directive 93/42/EEC, Regulation (EU) 2017/745 and ISO 13485 Medical Devices Quality Management Systems.

Classification: CE: Class I

FDA: Class II

The information contained within this manual must not be reproduced in whole or part without the manufacturer's prior written approval. As part of our policy for continued product development we the manufacturer reserve the right to make changes to specifications and other information contained in this document without prior notice.

This IFU is also available on the Keeler UK and Keeler USA websites.

Copyright © Keeler Limited 2021. Published in the UK 2021.

1. INDICATIONS FOR USE

Intended use / purpose of instrument

The Keeler Loupes and Frames are intended for use in various clinical settings including but not limited to users in dental practices, ophthalmology clinics, operating theatres, hospitals, veterinary practices, private consulting rooms and surgical environments and used by trained students and experienced professionals.

1.1 WARNINGS AND CAUTIONS

Please note that the proper and safe functioning of our instruments is only guaranteed if both the instruments and their accessories are exclusively from Keeler Ltd.

Observe the following precautions in order to ensure safe operation of the Loupe Frame.



WARNINGS

- Please note that the proper and safe functioning of our instruments is only guaranteed if both the instruments and their accessories are exclusively from Keeler Ltd .
- Check your Keeler product for signs of transport / storage damage prior to use.
- Do not use if the product is visibly damaged and periodically inspect it for signs of damage or misuse.
- US Federal Law restricts this device to sale by or on the order of a physician or practitioner.
- This device is intended to be used only by suitably trained and authorised healthcare professionals.
- This product should not be immersed in fluids.



CAUTION

- Use only genuine Keeler approved parts and accessories or device safety and performance may be compromised.
- The product has been designed to function safely when at an ambient temperature between +10°C and +35°C.
- Keep out of the reach of children.
- To prevent condensation from forming, allow instrument to come to room temperature before use.
- For indoor use only (protect from moisture).
- There are no user serviceable parts inside. Contact authorised service representative for further information.
- Failure to carry out recommended routine maintenance as per the instructions in this IFU may reduce the operational lifetime of the product.
- At product end of life dispose of in accordance with local environmental guidelines (WEEE).

1.2 CONTRAINDICATION

There is no restriction to patient population this device can be used with other than those outlined in the contraindications stated below.

The Keeler Loupe and Frame restrict the user's field and depth of view, forcing them to adopt awkward head-neck postures in order to see clearly.

2. CLEANING AND DISINFECTION INSTRUCTIONS

Only manual non-immersion cleaning as described should be used for this instrument. Do not autoclave or immerse in cleaning fluids.


1. Wipe the external surface with a clean absorbent, non-shedding cloth dampened with a de-ionised water / detergent solution (2% detergent by volume) or water / isopropyl alcohol solution (70% IPA by volume). Avoid optical surfaces.
2. Ensure that excess solution does not enter the instrument. Use caution to ensure cloth is not saturated with solution.
3. Surfaces must be carefully hand-dried using a clean non-shedding cloth.
4. Safely dispose of used cleaning materials.

3. SETTING UP AND USING THE LOUPE FRAME

1. Fit the Loupe unit to the frame via the pin mounting system as shown.



2. The protective shield is fitted by pushing the Shield into position as shown.

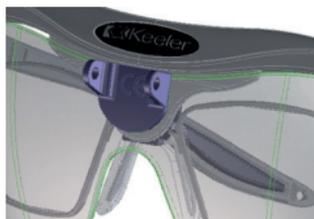

Protective shield in correct position

3. The Prescription Lens Frame is fitted if required as shown.



Protective shield in correct position

4. To adjust the nose bridge, turn the locking cam to a position as shown to release the Nose Bridge. The Nose Bridge can then be slid to the required position.



5. Lock the Nose Bridge in position by turning the locking cam to a position as shown.




6. To securely fit the frame to the user's head, slide the toggle indicated to a comfortable position



4. TECHNICAL SPECIFICATIONS

Environmental Conditions:

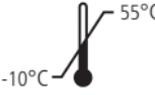USE	
  	
Shock (without packing)	10 g, duration 6 ms
STORAGE CONDITIONS	
  	
TRANSPORT CONDITIONS	
  	
Vibration, sinusoidal	10 Hz to 500 Hz: 0.5g
Shock	30 g, duration 6 ms
Bump	10 g, duration 6 ms

5. SPARE PARTS AND ACCESSORIES

Item	Part Number
Headband for surgical Loupes	2199-P-7259
Head strap	2199-P-7523
Hi-tech lens cloth in wallet	2199-P-7136

6. WARRANTY

No user serviceable parts – all preventative maintenance and servicing must only be performed by authorised Keeler representatives.

Your Keeler product is guaranteed for 3 years and will be replaced, or repaired free of charge subject to the following:

- Any fault due to faulty manufacture.
- The instrument and accessories have been used in compliance with these instructions.
- Proof of purchase accompanies any claim.

7. PACKAGING AND DISPOSAL INFORMATION

Disposal of old electrical and electronic equipment



This symbol on the product or on its packaging and instructions indicates that this product shall not be treated as household waste.

To reduce the environmental impact of WEEE (Waste Electrical Electronic Equipment) and minimise the volume of WEEE entering landfills we encourage at product end of life that this equipment is recycled and reused.

If you need more information on the collection reuse and recycling then please contact B2B Compliance on 01691 676124 (+44 1691 676124). (UK only).

Any serious incident that has occurred in relation to the device must be reported to the manufacturer and the competent authority of your Member State.

Contact



Manufacturer

Keeler Limited
Clewer Hill Road
Windsor
Berkshire
SL4 4AA UK

Freephone 0800 521251

Tel +44 (0) 1753 857177

Fax +44 (0) 1753 827145

USA Sales Office

Keeler USA
3222 Phoenixville Pike
Building #50
Malvern, PA 19355 USA
Toll Free 1 800 523 5620
Tel 1 610 353 4350
Fax 1 610 353 7814

India Office

Keeler India
Halma India Pvt. Ltd.
Plot No. A0147, Road No. 24
Wagle Industrial Estate
Thane West – 400604, Maharashtra
INDIA
Tel +91 22 4124 8001

China Office

Halma China Group
名称：沃迈（上海）机电有限公司
地址：上海市闵行区金都路1165弄
123号23幢一号厂房三层B座
电话：021-6151 9025



Visiometrics, S. L., Vinyals, 131
08221 Terrassa, Spain